

## Cook® Flexipet® Pipettes

Flexible polycarbonate pipettes used for denuding as well as gamete and embryo manipulation. Designed to resist breaking or scratching.

Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Description
<b>Denuding Pipettes—5 Vials of 10 (50 Total)</b>			
G26711	K-FPIP-1130-10BS-5	130	oocyte/embryo manipulation
G26712	K-FPIP-1140-10BS-5	140	oocyte/embryo manipulation
G26713	K-FPIP-1170-10BS-5	170	oocyte/embryo manipulation
<b>Manipulation Pipettes—1 Vial of 10 (10 Total)</b>			
G18878	K-FPIP-1300-10BS	300	blastocyst handling
<b>Manipulation Pipettes—5 Vials of 10 (50 Total)</b>			
G46020	K-FPIP-1300-10BS-5	300	blastocyst handling
G26057	K-FPIP-1600-10BS-5	600	manipulation of oocyte-cumulus complex



## Cook® Flexipet® Adjustable Handle


Used with Cook Flexipet pipettes for handling gametes and embryos.

- Aspiration volume can be preset from 0.25  $\mu$ L to 3.0  $\mu$ L.
- Provides the fine control necessary for oocyte, embryo and/or blastomere manipulation.

Global Product Number	Order Number	Description
G18674	K-MPH-1000	adjustable handle accepts all sizes of Flexipet pipettes


### Accessories

G18904	K-FRO-1000	one small O-ring, one large O-ring, spacer and collet
--------	------------	---

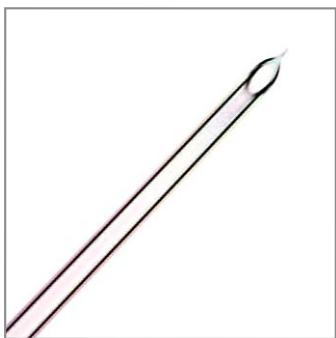

## Cook® Precision Micro-Injection Pipettes

Used for the intracytoplasmic single-sperm injection of oocytes.

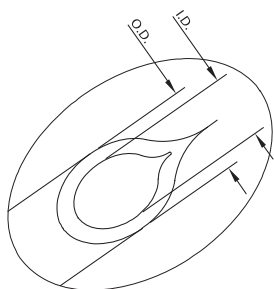


Pipettes specifically designed for ICSI, with parallel walls for precision sperm control and well-defined bevels to aid zona pellucida puncture. These pipettes are supplied in two main groups: the Standard "1" and the Precision "3" series. The Precision "3" series pipettes have parallel walls to better control sperm movements, while a stiff, spiked tip facilitates easy insertion through the zona pellucida. A longer bevel length provides more precise injection.

Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Outer Diameter $\mu\text{m}$	Distal Tip Angle	Unit of Sale
G32813	K-MPIP-1020	5.0	7	20°	box of 10
G32814	K-MPIP-1030	5.0	7	30°	box of 10
G18090	K-MPIP-1035	5.0	7	35°	box of 10
G32821	K-MPIP-3330	5.5	7	30°	box of 10
G26684	K-MPIP-3335	5.5	7	35°	box of 10




## Cook® Precision Holding Pipettes

Used to hold the oocyte in position by the application of vacuum when performing intracytoplasmic single-sperm injection or other micro-manipulation techniques.

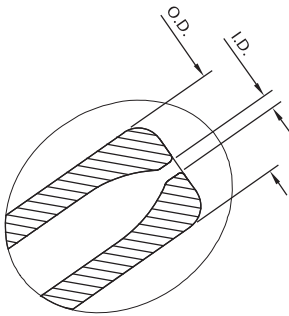


The pipettes are flame polished. The "2" series has a slightly larger inner and outer diameter to improve the stability of the oocyte during ICSI procedures. The "3" series of holding pipettes has a larger inner and outer diameter to provide greater stability and holding power to embryos during procedures such as assisted hatching and biopsy.

Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Outer Diameter $\mu\text{m}$	Distal Tip Angle	Unit of Sale
G32801	K-HPIP-1030	15	75	30°	box of 10
G32802	K-HPIP-1035	15	75	35°	box of 10
G32806	K-HPIP-2130	20	110	30°	box of 10
G32807	K-HPIP-2135	20	110	35°	box of 10
G32810	K-HPIP-3330	30	120	30°	box of 10
G32812	K-HPIP-3335	30	120	35°	box of 10

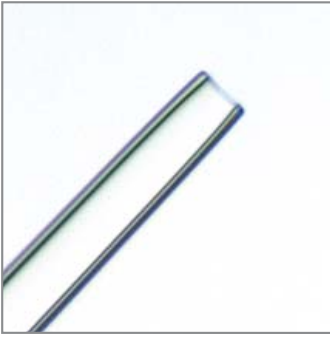

## Cook® Blastomere Aspiration Pipettes

Used for the aspiration of blastomeres to diagnose genetic disorders prior to embryo selection.

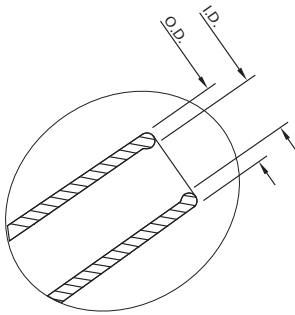


The pipette tip is cut straight and flame polished smooth. Parallel walls facilitate smooth aspiration.

Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Outer Diameter $\mu\text{m}$	Distal Tip Angle	Unit of Sale
G32795	K-EBPH-3035	30	42	35°	box of 10
G32799	K-EBPH-3535	35	49	35°	box of 10



## Cook® Polar Body Biopsy Pipette

Used for the aspiration of the polar body to diagnose genetic disorders prior to embryo selection.

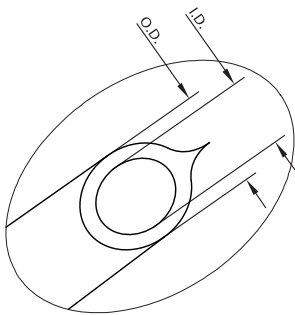


- The pipette tip is beveled with a small spike to aid zona pellucida insertion.
- The pipette has parallel walls to facilitate smooth insertion and aspiration of polar body.

Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Outer Diameter $\mu\text{m}$	Distal Tip Angle	Unit of Sale
G32822	K-PBBP-2035	20	28	35°	box of 10



## Cook® Zona Drilling Pipettes

Used to make a hole in the zona pellucida to enable assisted hatching and/or blastomere aspiration.

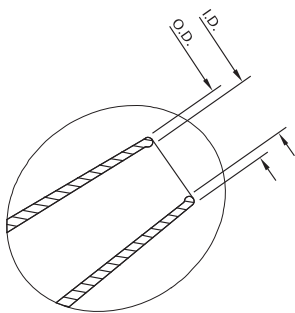

The pipette is designed to facilitate acid manipulation.


Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Outer Diameter $\mu\text{m}$	Distal Tip Angle	Unit of Sale
G32790	K-AHP-1035	5	10	35°	box of 10
G32794	K-AHP-2035	10	14	35°	box of 10



### Cook® Partial Zona Dissection Pipette

Used to create a slit or slits in the zona pellucida to enable assisted hatching and/or blastomere aspiration.

This pipette has a very fine point to create as smooth a cut as possible.

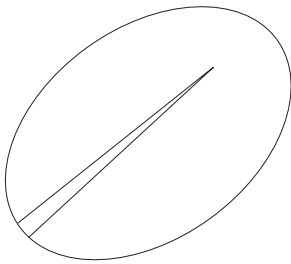


Global Product Number	Order Number	Distal Tip Angle	Unit of Sale
G18095	K-PZDP-1035	35°	box of 10



### Cook® Testicular Sperm Extraction Pipette

Used to extract sperm from the testicular tissue sample prior to intracytoplasmic sperm injection (ICSI).

The inner diameter enables easy extraction of sperm from the tissue sample. The parallel walls facilitate smooth aspiration. A Cook Precision Micro-Injection Pipette is then used for the sperm injection.

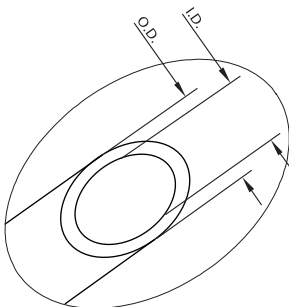

Global Product Number	Order Number	Inner Diameter $\mu\text{m}$	Outer Diameter $\mu\text{m}$	Distal Tip Angle	Unit of Sale
G32823	K-TSEP-1035	7	9	35°	box of 10