



IPC[®] System

INTEGRATED POWER CONSOLE



Innovation *that Transforms*

Inspired by Surgeons and Nurses in the Field

The Integrated Power Console is the
ONLY powered ENT system to offer:

- The widest range of application-specific products
- IntelliFlow® remote control for irrigation
- 0.5 ml/minute irrigation flow rates for working in the facial recess
- Integrated facial nerve monitoring with the Visao® Otologic Drill and Stim Bur Guard*
- High-speed, application-specific curved burs for otologic procedures
- A complete anterior skull base surgery solution

*For use with Medtronic's NIM ® 2.0 and 3.0 nerve monitoring systems

Integrated Power Console (IPC®) System *Innovation*

The new Integrated Power Console delivers innovation that transforms your O.R. Its design was inspired by feedback from you, the surgeons and nurses in the field. Our IPC® system is the only powered ENT surgery console with the most expansive range of application-specific products. From inferior turbinoplasties to complex skull base cases, the IPC platform is your ENT surgery system.

*Integrated facial nerve monitoring with the Visao® Otologic Drill and Stim Bur Guard**

Easy-to-use touchscreen menu for each handpiece

IntelliFlow® remote control for irrigation

Maximum visualization for anterior skull base

Connect multiple handpieces at once



*For use with Medtronic's NIM® 2.0 and 3.0 nerve monitoring systems

ion that Transforms

Complete Integration

Power all of your favorite Medtronic ENT handpieces

- Straightshot® M4 Microdebrider
- Endo-Scrub® 2 Lens Cleaning Sheaths
- Skeeter® Otologic Drill
- Visao® High-Speed Otologic Drill



Flexible Irrigation

On-drill irrigation reinvented

- Disposable IntelliFlow® remote allows precise flow control from sterile field
- 0.5 ml/minute irrigation flow rates for working in the facial recess
- Snap-on cassette for easier set-up
- Inner diameter of tubing enables continuous flow, even at low flow rates



Multifunction Footpedal

Ergonomic control at your feet


- Three buttons provide maximum control
- Built-in footrests and ergonomic angle for comfort
- Hoop locks firmly in place, then collapses easily for storage
- Non-slip design and compact size is ideal for your O.R.



Indigo™

High Speed Otologic Drill

- High speed and ample torque for fast debulking of the mastoid
- Available with straight and angled attachments
- Improved surgical access and visibility in narrow otologic spaces
- Works with a wide variety of Medtronic Xomed otology burs
- Grooves on attachment nose design allows for more stable irrigation

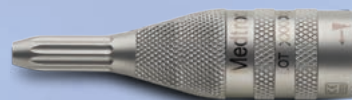

Anterior Skull Base Surgery

The complete solution for superior visualization


- Curved burs and blades allow better visibility and access
- Integrated suction keeps surgical site clear

The IPC platform is the ONLY powered ENT system that provides the widest range of application-specific products.

Stim Bur Guard for the Visao® Otologic Drill

Integrated facial nerve monitoring

- Provides 1-3 mm advance warning of facial nerve proximity
- Helps reduce the risk of patient injury for added peace-of-mind
- Works with our NIM® intraoperative nerve monitor*

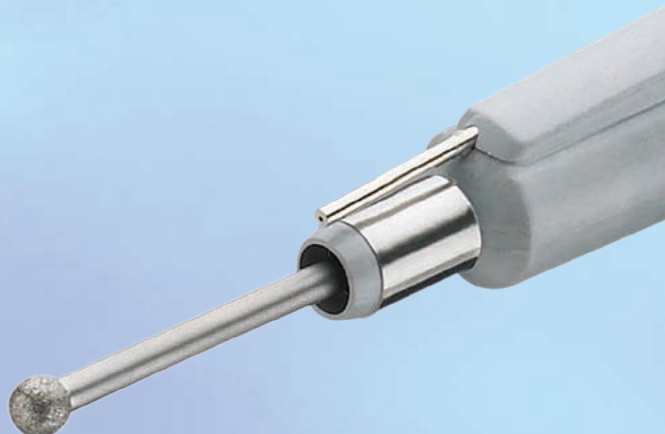

*For use with Medtronic's NIM® 2.0 and 3.0 nerve monitoring systems

Ordering Information

IPC® System and Endo-Scrub® 2 Accessories

Product	Description	Qty
1898001	Integrated Power Console	1
1898430	Integrated Power Console Multifunction Footpedal	1
1897510	Basket	1
1852000	Footswitch Endo-Scrub 2	1
1991015	Finger Switch Endo-Scrub 2	1

IPC® System Power Cords

Product	Description	Qty
1895820	Power Cord, Standard, North America, 3M	1
1895821	Power Cord, Standard, North America, 6M	1
1895824	Power Cord, UK/Ireland, 6M	1
1895825	Power Cord, Continental Europe, 6M	1
1895823	Power Cord, Japan	1

IPC® System Manuals

Product	Description	Qty
1898851	IPC Manual (English Only)	1
1898851A	IPC Manual (EL, EN, ES, FR, PT)	1
1898851B	IPC Manual (DE, EN, FR, IT, NL)	1
1898851C	IPC Manual (DA, EN, FI, NO, SV)	1
1898851D	IPC Manual (CS, EN, HU, PL, TR)	1
1898851E	IPC Manual (EN, RU) - CD only	1

Straightshot® M4 Microdebrider Components

Product	Description	Qty
1898200T	XPS® Straightshot® M4 High-Speed Microdebrider	1
1898400	XPS Straightshot M4 Instrument Tray	1

For further information, please call Medtronic ENT at 800-874-5797 or 904-296-9600. You may also consult our website at www.MedtronicENT.com.

Medtronic ENT

Medtronic USA, Inc.
6743 Southpoint Drive North
Jacksonville, FL 32216 USA
www.MedtronicENT.com
Telephone: (904) 296-9600
Toll free: (800) 874-5797
Fax: (904) 296-9666

International Telephone Numbers

Australia 1-800-668-670	Korea 82-2-3404-3600
Belgium 32-2456-09-00	Lebanon 961-1-370-670
Canada 1-800-217-1617	Luxembourg 32-2456-09-09
China 86-21-50800998	Netherlands 31-45-566-8000
France 33-470-679-800	Poland 48-22-465-6942
Germany 49-211-5293-209	Singapore 65-6776-6255
Hong Kong 852-2919-1312	Spain 34-91-625-05-40
India 91-22-26836733	UK 44-1923-212-213
Italy 39-02-24137-324	USA 1-904-296-9600
Japan 81-6-4795-1506	

Visao® High-Speed Otologic Drill Components

Product	Description	Qty
3334800	Visao High-Speed Otologic Drill	1
3334625	Visao High-Speed Otologic Drill Extended Bur Guard without Irrigation	1
3334635	Visao High-Speed Otologic Drill Extended Bur Guard with Irrigation	1
3334895	Visao High-Speed Otologic Drill Instrument Tray	1
3155617	Bur Brush	1
3112500	Micro-Craft Chuck Brush Cleaner	60
3318603	Visao Irrigation Tubing with IntelliFlow Remote Control	1
3318605	Suction/Irrigation Connector Set	5
3318606	Visao Coolant Tubing	2
3318601	Stim Bur Guard with Irrigation	1
3318602	Stim Bur Guard without Irrigation	1

Skeeter® Components

Product	Description	Qty
3055601	Skeeter Ultralite Handpiece	1
3155615	Skeeter Instrument Tray	1
3052460	Skeeter Otologic Metric Scale	1
3155641	Skeeter Oto-Flex Bur, Orange/White, Diamond, 1 mm	1
3155631	Skeeter Oto-Flex Bur, Orange, Carbide, 1 mm	1
3155647	Skeeter Oto-Flex Bur, Green/White, Diamond, 0.7 mm	1
3155632	Skeeter Oto-Flex Bur, Red, Carbide, 2.3 mm	1
3155648	Skeeter Oto-Flex Bur, Yellow/White, Diamond, 0.8 mm	1
3155642	Skeeter Oto-Flex Bur, Red/White, Diamond, 2.3 mm	1
3155646	Skeeter Oto-Flex Bur, Blue/White, Diamond, 0.6 mm	1
3155643	Skeeter Oto-Flex Bur, Gray/White, Diamond, 1.4 mm	1
3155645	Skeeter Oto-Flex Bur, Brown/White, Diamond, 1.8 mm	1
3155635	Skeeter Oto-Flex Bur, Violet, Carbide, 0.5 mm	1

Indigo™ Components

Product	Description	Qty
1845000	Indigo High-Speed Otologic Drill	1
1845010	Indigo Otologic Straight Attachment	1
1845020	Indigo Otologic Angled Attachment	1
1845040	Indigo Instrument Tray	1
1845030	Indigo Tubing	

This literature is intended for the exclusive use of physicians. Rx only.
©2013, Medtronic, Inc. *Registered mark of Medtronic, Inc.
UC201400184 EN 06.13